FLOCARE® BENCHMARK® NI TUBE

Manufacturer: Nutricia Medical Devices BV, Taurusavenue 167, 2132 LS Hoofddorp, The Netherlands, www.nutriciaflocares.com
Production location: Nutricia Pharmaceutical Wuxi Co., Ltd., P.R. China

STERILE

EO

MD

200344_D
May 2022

FLOCARE® BENCHMARK® NI TUBE

EN • INSTRUCTIONS FOR USE

BENCHMARK® NI TUBE

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BENCHMARK® TRANSNASALE DÜNNDARMSONDE

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BENCHMARK® NASO-INTESTINALE SONDE

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SONDE NASO-INTESTINALE BENCHMARK®

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SONDA NASOINTESTINAL BENCHMARK®

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SONDINO NASOINTESTINALE BENCHMARK®

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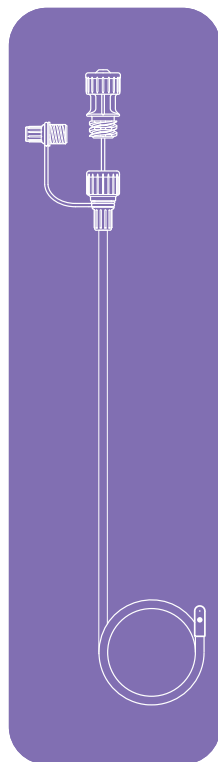
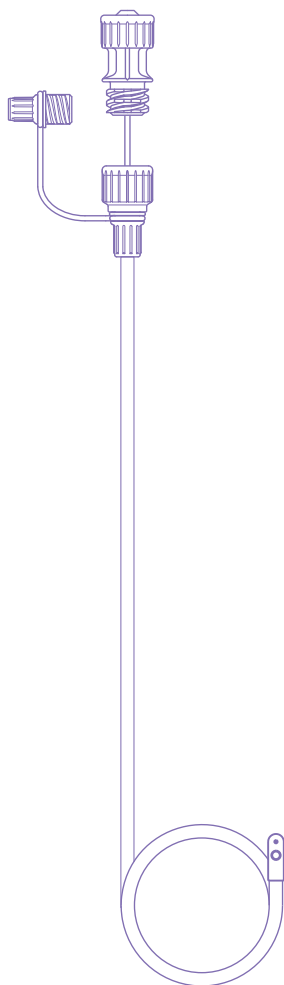
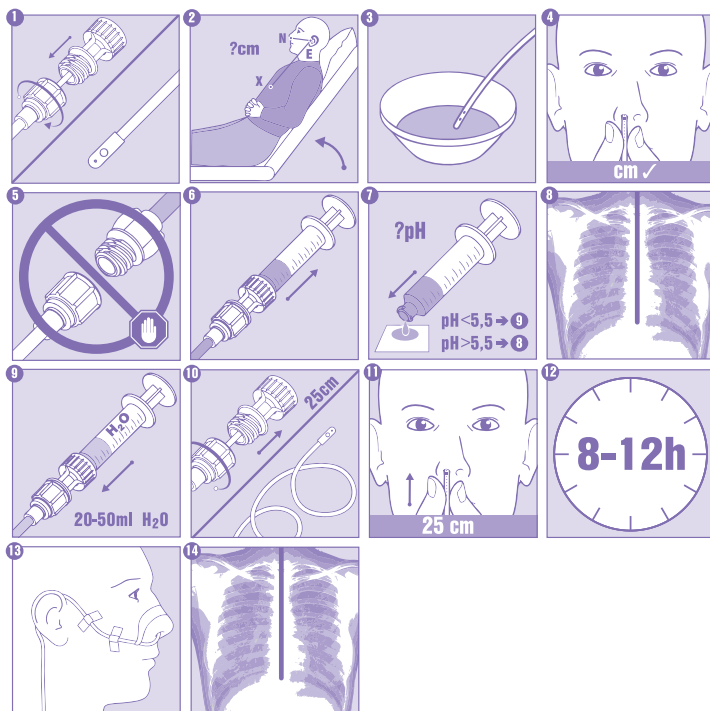
EN • BENCHMARK® NI TUBE

PRODUCT DESCRIPTION

The Flocare Benchmark® Naso-Intestinal feeding tube is a transnasal intestinal tube with a Benchmark® spiral tip designed for the administration of enteral feeding.

The tube has the following characteristics:

Tube material: white full opaque polyurethane tubing

External diameter: CH8 or CH10

Length: 145 cm

Tip: bullet tip with 2 lateral outlets

Hydromer® coating on the tip

Benchmark® spiral (patented): 2.5 loops with a diameter of approx. 3cm and a length of approx. 23 cm. These loops are straightened with a guidewire during insertion. After removal of the guidewire, the spiral eases transpyloric passage and allows optimal fixation in the small intestine.

Guidewire: A prelubricated metallic guidewire, partially introduced into the tube.

Connector: ENFit connector.

INDICATIONS FOR TUBE FEEDING OR ENTERAL NUTRITION

Enteral nutrition or tube feeding is indicated to patients at nutritional risk or malnourished who:

- Cannot meet their nutrient requirements by normal dietary intake and;
- Have a functional gastrointestinal tract

INDICATIONS FOR NASO-INTESTINAL FEEDING TUBES USE

For short-term administration use – up to 6 weeks.

Indication for therapeutic reasons:

- Enteral feeding with functional intestinal tract with impaired gastric emptying
- Post-pyloric feeding directly into the small intestine
- Patients with an increased risk of enteral diet regurgitation and/or pulmonary aspiration.

It is recommended to always deliver feed in the small intestine using an enteral feeding pump to control the flow rate.

CONTRA-INDICATIONS

Absolute contraindications:

- Gastrointestinal obstruction
- Perforation/ leakage

Relative contraindications:

- Severe decreased gut function i.e.: severe malabsorption, paralytic or mechanical ileus
- Motility disorders
- Acute abdomen
- Severe Crohn's disease

Notice: any serious incident that has occurred in relation to this device is to be reported to Nutricia and the local competent authority.

TUBE PLACEMENT

Naso-intestinal feeding tubes should be inserted by trained healthcare professionals because of the risk of misplacement and esophageal or pulmonary perforation.

Instructions for use for patients with (limited) gastric motility:

1. Advise the patient of the placement procedure. Fully insert the guide wire in the tube and ensure that it is firmly attached to the connector (fig.1).

2. Place the patient in a comfortable sitting or semi-recumbent position. To determine the length of Bengmark® tube to be inserted, measure the distance from the tip of the Nose (N), to the Ear (E), and then to the Xiphoid process of the sternum (X) (fig.2). Mark the tube at this point. Make 2 additional marks, at 25 cm and at 50 cm after the first mark.

3. Lubricate the tip and first cm of the tube with water, to facilitate insertion (fig.3).

4. Inspect the patient's nostrils for a deviated septum or visible nasal polyps; use the larger nostril and pass the tube along the bottom of the nostril (fig.4). When the tube reaches the throat, bend the patients' head forward and ask them to swallow as often as possible. Push the tube gently, never force it. Avoid inadvertent tube placement in the trachea! Continue until the first pre-measured mark is reached. Giving something to drink or to eat (if possible) can stimulate the migration of the tube.

Warning: Do not connect the Bengmark® NI tube with the feeding set at this stage (fig 5). Before correct placement is confirmed do not administer any fluids through the tube:

- Do not flush the tube

- Do not lubricate the internal guidewire

5. Confirm position of the tube:

(fig 6) Gastric content aspiration

Obtain an aspirate of the gastric contents through the ENFit connector on the guidewire handle.

(fig 7) pH measurement

Measure the pH value of the aspirate. A pH $\leq$ 5.5 indicates a gastric placement. Local guidelines can differ in terms of the acceptable pH range for confirming NG tube placement. Always consult your local guidelines

(fig 8) X-ray confirmation

In case no aspirate is obtained, or the measured pH value is 5.5 or higher; Confirm gastric placement by X-ray.

Warning: Please follow the local radiation protection procedures.

6. After checking the position of the tube (fig. 6, 7, 8) inject at least 20-50 ml H₂O to remove gastric content from the tube. (fig.9).

7. Pull $\pm$ 25cm of guidewire out of the tube (fig.10). Then continue passing the tube until you reach the second mark (fig.11). Now pull the guidewire out of the tube completely.

8. Do not yet secure the tube to the nose, use the ear lobe as a temporary fixation point.

9. Within 8-12 hours gastric motility will have moved the tube through the pylorus (fig.12). Secure the tube once the third mark on the tube is at the patients' nose (fig. 13).

10. Before administrating any enteral nutrition, check if the tube is positioned correctly by X-ray (fig.14).

Warnings:

- Coughing may indicate passage of the tube into the trachea. If this is suspected, remove the tube and reinsert it once the patient is comfortable.

- If resistance is encountered, remove the tube. Notify physician.

- Particular care should be taken if an endotracheal tube is in place, as it may tend to guide the feeding tube into the trachea.

- Failure of transpyloric passage: following appropriate consultation with a clinician, the use of gut motility agents may be considered and should be given to the patients 10 min prior to passing the tube.

TUBE PLACEMENT USING FLUOROSCOPY, ENDOSCOPY OR 2-STAGE X-RAY

For patients who are unconscious, intubated, sedated, ventilated, have minimal or absent gag reflex and/or are otherwise debilitated or have no stomach motility at all, one of the following procedures may be considered as an alternative tube placement.

1. Tube Placement using Fluoroscopy

2. Tube placement using Endoscopy

Use standard fluoroscopy/endoscopy procedures for feeding tubes intragastric placement.

Additional steps for placing the NI tube to pass through the pylorus:

Grasp the tube with forceps of endoscope, and advance endoscope and tube through pylorus

Place the tube as far as possible into the small intestine, preferably beyond ligament of Treitz

Leave the tube with forceps in place and withdraw endoscope as far as possible

Flush the tube with water via guidewire connector to stiffen tube

Remove guidewire and withdraw forceps and endoscope carefully

Verify tube placement with x-ray

3. Tube Placement using 2-stage X-Ray

Stage 1 : confirm the tube has not entered the bronchial tree

Stage 2 : confirm final placement

MEDICATION ADMINISTRATION AND AFTER CARE

Caregivers should apply a risk-based approach before the nasoenteric route is selected to administer medication to the patient. Appropriate ancillaries including syringes should be used for medication administration through enteral feeding tubes using ENFit™ connectors of a recognized standard (ISO80369-3) in order to avoid misconnection errors.

Check the position of the tube (cm mark at the nostril remains the same):

- Each time the tube feed is changed.
- When in doubt whether the tube is in the correct position.

Check the fixation tape regularly.

GENERAL FLUSHING AND SYRINGE USE POLICY

- Routinely flush the tube with water before and after the administration of feed or medication and at least every 8 hours to prevent blockage of the tube.

- As practicable, use an enteral syringe (minimum size of 20mL) for flushing the tube and always flush in a controlled manner.

- For medication administration directly into the tube, always administer using an appropriate size enteral syringe for the volume of medication but ensure delivery is in a controlled manner. In case more than one medication is administered at one time; make sure the tube is flushed thoroughly between types.

- Never use excessive force. If there is any resistance, stop the procedure and consult your clinician.

REMOVAL

Flocare Benchmark® Naso-intestinal feeding tube is recommended for short-term tube feeding up to 6 weeks. Replacement or removal of the product may differ from this recommendation, based on medical scrutiny, viability of the product and/or potential patient inconveniences and is to be executed under medical scrutiny and supervision. Tube replacement should be accomplished in case of tube breakage, occlusion, dislodgement or degradation.

The Benchmark® Naso-intestinal Tube can be removed by gently pulling the tube out of the patient's nose. Removal should always be performed under medical supervision.

After use, dispose of product and packaging in accordance with hospital, administrative and/or local policy.

POTENTIAL COMPLICATIONS

The commonly reported complications of nasoenteric feeding tubes are respiratory due to tube misplacement. Other potential mechanical complications are listed below:

- incorrect intubation,
- tube obstruction,
- misconnections,
- tube breakage/kinking,
- tube dislocation,
- irritation/pain of oesophagus, pharynx or around nostrils of the nose,
- pressure injury related to fixation.

General potential complications for enteral feeding are gastro-intestinal, metabolic and infectious.

SINGLE PATIENT USE

Do not reuse this product - Do not reuse, reprocess or re-sterilize this device. Reuse, reprocessing or re-sterilization may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death. Reuse, reprocessing or re-sterilization may also create a risk of contamination of the device and/or cause patient infection or cross-infection, including, but not limited to, the transmission of infectious disease(s) from one patient to another. Contamination of the device may lead to injury, illness or death of the patient.

GBR: Nutricia Ltd, White Horse Business Park, TROWBRIDGE, Wiltshire, BA14 0XQ, England. Tel: 01225 711 688, www.nutricia.co.uk Patient/Carer Information Number (UK only) Tel: 08457 623 653.

IRL: Nutricia Ireland Ltd., Block 1, Deansgrange Business Park, Deansgrange, Co. DUBLIN, Ireland. www.nutricia.ie

AUS: Nutricia Australia Pty Ltd, Lvl 4, Bld D, 12-24 Talavera Rd, Macquarie Park NSW 2113, Australia, 1800 060 051, www.nutriciamedical.com.au

NZL: Danone Nutricia NZ Ltd, 56-58 Aintree Avenue, Airport Oaks, Auckland 2022, New Zealand, 0800 636 228, www.nutriciamedical.co.nz.

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